


Performs best under pressure

Safety. Convenience. Reliability. Cost. These are the most important tourniquet issues that concern physicians, OR personnel, and administrators. With a special focus on these issues, Zimmer incorporates the most practical elements of today's technology to develop the A.T.S.[®] 750 Automatic Tourniquet System.

The A.T.S. 750 Tourniquet uses ambient air so there is no need for nitrogen tanks or wall hook-ups. Microprocessor controls provide accuracy and speed, allowing the system to react quickly to manual input changes or self-monitored conditions that require automatic adjustments.




Disposable and reusable cuffs come in a variety of sizes and styles designed to fit the largest and smallest of limbs.

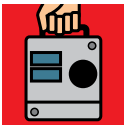

Small and Lightweight- the unit is compact, weighing in at only 6 pounds.



Audible/Visual Alarm- audible and visual alarms tell you when there is an unsafe condition.

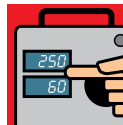

Self Calibration- the unit performs a self-calibration check and diagnostic test every time it is turned on.



Cuff Inflation Security- lockout feature ensures that the cuff stays inflated, even if the unit is turned off.



Touchpad switches for time and pressure controls are easy to set, and large LED displays are easy to read.




Hoses are easy to connect and use a positive locking mechanism to help prevent unintentional detachment.



The unit is compact, lightweight, and easy to use. For added reliability, a battery back-up helps ensure continuous power. And for added peace of mind, a visual and audible alarm is activated when an unsafe condition is detected, from a pressure change to a leaking cuff.

Put the pressure where it belongs


The A.T.S. 750 Tourniquet is a full-featured, automatic tourniquet system that is sensitive to the cost issues that all health care professionals face today. It is the perfect balance of features, technology, and practicality. By allowing you to control every aspect of the tourniquet without diverting concentration from other important intraoperative matters, the A.T.S. 750 Tourniquet helps take the pressure off you, and put it right where it belongs.

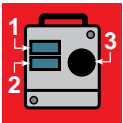

Simple Controls- controls are easy to learn and use.



Positive Locking Connectors- positive locking connectors help prevent unintentional detachment.

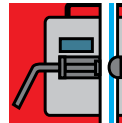

Mounting Bracket and Footpads- unit can be mounted on an I.V. pole or placed on a table top.



Rotary dial allows fast, accurate adjustments of time or pressure values. It can be used to change both current and default settings.

ORDERING INFORMATION

A.T.S. 750 TOURNIQUET SYSTEM

SYSTEMS

Prod. No.	Description
60-0750-101	A.T.S. 750 with Hose Only. No Stand, Basket, or Cuffs.
60-0750-102	A.T.S. 750 with Hose and A.T.S. Cylindrical Cuffs (60-7600-01/07)
60-0750-103	A.T.S. 750 with Hose, A.T.S. Cylindrical Cuffs (60-7600-01/07), Stand, and Basket
60-0750-104	A.T.S. 750 with Hose, Stand, and Basket. No Cuffs.

STERILE DISPOSABLE CUFFS

SINGLE PORT, SINGLE BLADDER (BOX 10)

PLC	Luer-Lock	PLC	Luer-Lock	Size
Without Sleeves	Without Sleeve	With Sleeve	With Sleeve	
60-7075-101	60-7075-01	60-7075-151	60-7075-51	8 in.
60-7075-102	60-7075-02	60-7075-152	60-7075-52	12 in.
60-7075-103	60-7075-03	60-7075-153	60-7075-53	18 in.
60-7075-104	60-7075-04	60-7075-154	60-7075-54	24 in.
60-7075-105	60-7075-05	60-7075-155	60-7075-55	30 in.
60-7075-106	60-7075-06	60-7075-156	60-7075-56	34 in.
60-7075-107	60-7075-07	60-7075-157	60-7075-57	42 in.

CUFFS WITH SLEEVES

SINGLE PORT, SINGLE BLADDER (BOX 1)

PLC	Luer-Lock	Size
With Sleeve	With Sleeve	
60-7075-153-01	60-7075-53-01	18 in.
60-7075-156-01	60-7075-56-01	30 in.

SINGLE PORT, DUAL BLADDER (BOX 10)

PLC	Luer-Lock	PLC	Luer-Lock	Size
Without Sleeves	Without Sleeve	With Sleeve	With Sleeve	
60-7085-101	60-7085-01	60-7085-151	60-7085-51	12 in.
60-7085-102	60-7085-02	60-7085-152	60-7085-52	18 in.
60-7085-103	60-7085-03	60-7085-153	60-7085-53	24 in.

CUFF SLEEVES ONLY (BOX 10)

Prod. No.	Description
60-8000-08	Sleeve for 8 in. Disposable Cuff
60-8000-12	Sleeve for 12 in. Disposable Cuff
60-8000-18	Sleeve for 18 in. Disposable Cuff
60-8000-24	Sleeve for 24 in. Disposable Cuff
60-8000-30	Sleeve for 30 in. Disposable Cuff
60-8000-34	Sleeve for 34 in. Disposable Cuff
60-8000-42	Sleeve for 42 in. Disposable Cuff

ACCESSORIES

Prod. No.	Description
60-2790-01	Cuff Hose
60-4022-01	Tourniquet Stand (Includes Base and Pole)
60-1908-01	Accessory Basket



REUSABLE CUFFS (BOX 1)

A.T.S. CYLINDRICAL CUFFS

SINGLE PORT, SINGLE BLADDER

Prod. No.	Size
60-7600-01	8 in.
60-7600-02	12 in.
60-7600-03	18 in.
60-7600-04	24 in.
60-7600-05	30 in.
60-7600-06	34 in.
60-7600-07	42 in.

A.T.S. CONTOUR CUFF

SINGLE PORT, SINGLE BLADDER

Prod. No.	Description	Size
60-7600-08	Large Contour, Thigh	42 in.

A.T.S. CYLINDRICAL CUFF

SINGLE PORT, DUAL BLADDER

Prod. No.	Size
60-7666-01	12 in.
60-7666-02	18 in.
60-7666-03	24 in.

Please refer to package insert for complete product information, including contraindications, warnings, precautions, and adverse effects.

Contact your Zimmer Representative or visit us at www.zimmer.com.



zimmer
Confidence in your hands™



A.T.S.[®] 750 AUTOMATIC TOURNIQUET SYSTEM

Safe, Convenient, Reliable


A.T.S.[®]
750
AUTOMATIC TOURNIQUET SYSTEM

mmHg



250

50

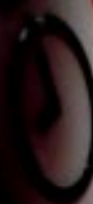

1/0